



ICO240

Safety at all times: High-end functions for the protection of cell cultures, bacteria cultures or tissue cultures.



On this page, you can find all the essential technical data on the Memmert CO₂ incubator ICO. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at sales@memmert.com.

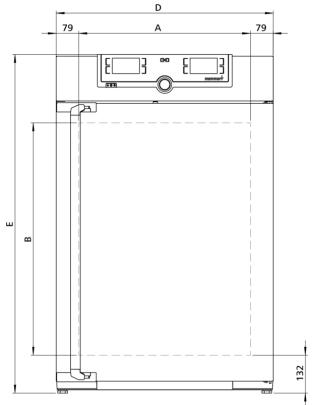

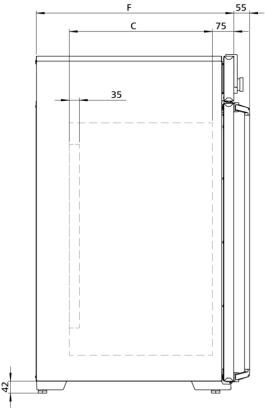

Temp	eratı	ıre
------	-------	-----

Working-temperature range	5 °C above ambient temperature up to +50 °C Standard sterilisation programme: 60 minutes at 180°C (without removing the sensors)
Setting temperature range	+18 to +50 °C
Setting accuracy temperature	0.1 °C
Temperature sensor	2 Pt100 sensors DIN Class A in 4-wire-circuit for mutual monitoring, taking over functions in case of an error
Temperature variation in chamber	at + 37 °C +/- 0.3 K
Temperature fluctuations with time	at 37 °C +/- 0.1 K

Humidity

Humidity control (standard)	Humidity limitation thanks to a Peltier element; when water dish is full and inserted, the Peltier element limits the value of relative humidity in the interior to 93 $\%$ rh +/- 2.5 $\%$	
Setting accuracy humidity	0.5 % rh	
Setting range active humidity control (with option K7)	40 to 97 % rh and rh-Off	

Control of standard components

CO ₂ control	Digital electronic CO ₂ control with dual beam NDIR system, with auto-diagnostic system and acoustic fault indication, barometric pressure compensation
Adjustment range CO ₂	0 to 20 % CO ₂
Setting accuracy CO ₂	0,1%
Variation in time CO ₂	+/- 0.2 % CO ₂
Adjustment range O ₂	1 to 20 % O ₂
Setting accuracy O ₂	0.1 % O ₂

Control technology

ControlCOCKPIT TwinDISPLAY. Adaptive multifunctional digital PID-microprocessor controller with 2 high TFT-colour displays.	
Language setting	German, English, Spanish, French, Polish, Czech, Hungarian, Italian, Chinese
Function SetpointWAIT	the process time does not start until the set temperature is reached
Adjustable parameters temperature (Celsius or Fahrenheit), CO ₂ , programme time, time zones, summertime/w	

_					
Col	mm	un	ıca	ıtion	۱

	EU LANGUED	
Interface	e Ethernet LAN, USB	
Documentation	programme stored in case of power failure	
Programming	AtmoCONTROL software on a USB stick for programming, managing and transferring programmes via Ethernet interface or USB port	
Safety		
AutoSAFETY	additionally integrated over- and undertemperature protection "ASF", automatically following the setpoint value at a preset tolerance range, alarm in case of over- or undertemperature, heating function is switched off in case of overtemperature, cooling function in case of undertemperature	
Autodiagnostic system	integral fault diagnostics for temperature and CO ₂	
Alarm	visual and acoustic	
Heating concept		
6 sides	large-area multi-function heating system on four sides with additional door and back heating to avoi condensation	

Standard equipment

Works calibration certificate	Works calibration certificate at +37 $^{\circ}$ C, 5 $^{\circ}$ CO $_2$ and 90 $^{\circ}$ rh (requires option K7); standard equipment for units with active humidity control	
Works calibration certificate	Works calibration certificate at +37 °C, 5 % $\rm CO_2$, 90 % rh and 10 % $\rm O_2$ (requires option K7 and option T6); standard equipment for units with $\rm O_2$ control	
Works calibration certificate	incl. works calibration certificate (measuring point chamber centre) at +37°C, 5 % $\rm CO_2$ for standard units	
Door	fully insulated stainless steel door with 2-point locking (compression door lock)	
Door	inner glass door with opening (8 mm Ø) to take gas sample	
Standard equipment	CO ₂ connection set: hose with coupling and clamp	
Standard accessories	Membrane filter (in order to remove impurities and pollutants, all incoming gases pass through a membrane filter before they reach the chamber)	
Internals	2 perforated stainless steel shelf/shelves	
Internals	1 stainless steel water dish	

Stainless steel interior

Dimensions	w _(A) x h _(B) x d _(C) : 600 x 810 x 500 mm (d less 35 mm for fan)
Interior	material 1.4301 (ASTM 304), corrosion resistant
Volume	241
Max. number of internals	12
Max. loading of chamber	140 kg
Max. loading per internal	15 kg

Textured stainless steel casing

Dimensions	w _(D) x h _(E) x d _(F) : 759 x 1180 x 691 mm (d +56mm door handle)
Housing	rear zinc-plated steel

Electrical data

Voltage	230 V, 50/60 Hz
Electrical load	approx. 1650 W
Voltage	115 V, 50/60 Hz
Electrical load	approx. 1650 W

Ambient conditions

Set Up	The distance between the wall and the rear of the appliance must be at least 15 cm. The clearance from the ceiling must not be less than 20 cm and the side clearance from walls or nearby appliances must not be less than 5 cm.	
Ambient temperature	10 °C to 35 °C	
Humidity rh	max. 70 %, non-condensing	
Altitude of installation	max. 2,000 m above sea level	
Overvoltage category	II	
Pollution degree	2	

Packing/shipping data

Transport information	The appliances must be transported upright
Customs tariff number	8419 8998
Country of origin	Federal Republic of Germany
WEEE-RegNo.	DE 66812464
Dimensions approx incl. carton	w x h x d: 840 x 1360 x 900 mm
Net weight	approx. 110 kg
Gross weight carton	approx. 145 kg

Standard units are safety-approved and bear the test marks





