


### **ICO50**

Safety at all times: High-end functions for the protection of cell cultures, bacteria cultures or tissue cultures.



On this page, you can find all the essential technical data on the Memmert CO<sub>2</sub> incubator ICO. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at <a href="mailto:sales@memmert.com">sales@memmert.com</a>.

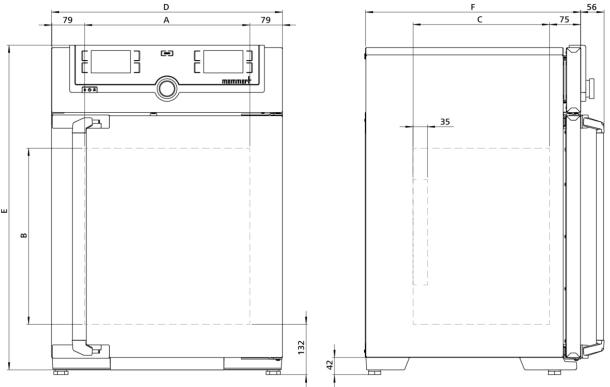

Te	m	pe	ra	tui	re
	•••	~	. ~	••••	_

Working-temperature range	5 °C above ambient temperature up to +50 °C
	Standard sterilisation programme: 60 minutes at 180°C (without removing the sensors)
Setting temperature range	+18 to +50 °C
Setting accuracy temperature	0.1 °C
Temperature sensor	2 Pt100 sensors DIN Class A in 4-wire-circuit for mutual monitoring, taking over functions in case of an error
Temperature variation in chamber	at + 37 °C +/- 0.3 K
Temperature fluctuations with time	at 37 °C +/- 0.1 K

# Humidity

Humidity control (standard)	Humidity limitation thanks to a Peltier element; when water dish is full and inserted, the Peltier element limits the value of relative humidity in the interior to 93 $\%$ rh +/- 2.5 $\%$
Setting accuracy humidity	0.5 % rh
Setting range active humidity control (with option K7)	40 to 97 % rh and rh-Off

# **Control of standard components**

CO <sub>2</sub> control	Digital electronic CO <sub>2</sub> control with dual beam NDIR system, with auto-diagnostic system and acoustic fault indication, barometric pressure compensation
Adjustment range CO <sub>2</sub>	0 to 20 % CO <sub>2</sub>
Setting accuracy CO <sub>2</sub>	0,1%
Variation in time CO <sub>2</sub>	+/- $0.2~\%~\mathrm{CO}_2$
Adjustment range O <sub>2</sub>	1 to 20 % O <sub>2</sub>
Setting accuracy O <sub>2</sub>	0.1 % O <sub>2</sub>

# **Control technology**

ControlCOCKPIT	TwinDISPLAY. Adaptive multifunctional digital PID-microprocessor controller with 2 high-definition TFT-colour displays.
Language setting	German, English, Spanish, French, Polish, Czech, Hungarian, Italian, Chinese
Function SetpointWAIT	the process time does not start until the set temperature is reached
Adjustable parameters	temperature (Celsius or Fahrenheit), ${\rm CO}_2$ , programme time, time zones, summertime/wintertime

_								
"	$\hat{}$	m	m	111	١ic	• • •	10	n
v	v			uı	IIL	٠aι	ıv	ш

Volume

Max. number of internals

Max. loading of chamber

Max. loading per internal

56 I

75 kg

15 kg

5

Interface	Ethernet LAN, USB
Documentation	programme stored in case of power failure
Programming	AtmoCONTROL software on a USB stick for programming, managing and transferring programmes via Ethernet interface or USB port
Safety	
AutoSAFETY	additionally integrated over- and undertemperature protection "ASF", automatically following the setpoint value at a preset tolerance range, alarm in case of over- or undertemperature, heating function is switched off in case of overtemperature, cooling function in case of undertemperature
Autodiagnostic system	integral fault diagnostics for temperature and ${\rm CO}_2$
Alarm	visual and acoustic
Heating concept	
6 sides	large-area multi-function heating system on four sides with additional door and back heating to avoid condensation
Standard equipment	
Works calibration certificate	Works calibration certificate at +37 $^{\circ}$ C, 5 $^{\circ}$ CO $_2$ and 90 $^{\circ}$ rh (requires option K7); standard equipment for units with active humidity control
Works calibration certificate	Works calibration certificate at +37 °C, 5 % $\rm CO_2$ , 90 % rh and 10 % $\rm O_2$ (requires option K7 and option T6); standard equipment for units with $\rm O_2$ control
Works calibration certificate	incl. works calibration certificate (measuring point chamber centre) at +37°C, 5 $\%~\mathrm{CO_2}$ for standard units
Door	fully insulated stainless steel door with 2-point locking (compression door lock)
Door	inner glass door with opening (8 mm $\varnothing$ ) to take gas sample
Standard equipment	CO <sub>2</sub> connection set: hose with coupling and clamp
Standard accessories	Membrane filter (in order to remove impurities and pollutants, all incoming gases pass through a membrane filter before they reach the chamber)
Internals	1 perforated stainless steel shelf/shelves
Internals	1 stainless steel water dish
Stainless steel interior	
Dimensions	w <sub>(A)</sub> x h <sub>(B)</sub> x d <sub>(C)</sub> : 400 x 425 x 330 mm (d less 35 mm for fan)
Interior	material 1.4301 (ASTM 304), corrosion resistant

### Textured stainless steel casing

Dimensions	w <sub>(D)</sub> x h <sub>(E)</sub> x d <sub>(F)</sub> : 559 x 795 x 521 mm (d +56mm door handle)
Housina	rear zinc-plated steel

#### **Electrical data**

Voltage	230 V, 50/60 Hz
Electrical load	approx. 1100 W
Voltage	115 V, 50/60 Hz
Electrical load	approx. 1100 W

#### **Ambient conditions**

Set Up	The distance between the wall and the rear of the appliance must be at least 15 cm. The clearance from the ceiling must not be less than 20 cm and the side clearance from walls or nearby appliances must not be less than 5 cm.
Ambient temperature	10 °C to 35 °C
Humidity rh	max. 70 %, non-condensing
Altitude of installation	max. 2,000 m above sea level
Overvoltage category	II .
Pollution degree	2

### Packing/shipping data

Transport information	The appliances must be transported upright
Customs tariff number	8419 8998
Country of origin	Federal Republic of Germany
WEEE-RegNo.	DE 66812464
Dimensions approx incl. carton	w x h x d: 730 x 950 x 640 mm
Net weight	approx. 55 kg
Gross weight carton	approx. 74 kg

#### Standard units are safety-approved and bear the test marks





