



UF30m

The heating oven Um is a Class I medical device.

This oven made of high-quality hygienic, easy to clean stainless steel leaves nothing to be desired in terms of ventilation technology, control technology, overtemperature protection and perfectly tuned heating technology.

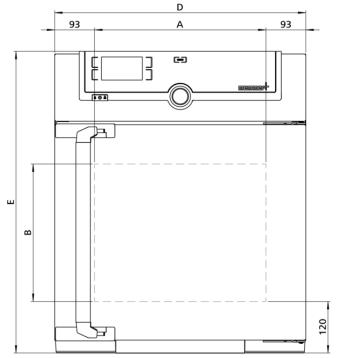

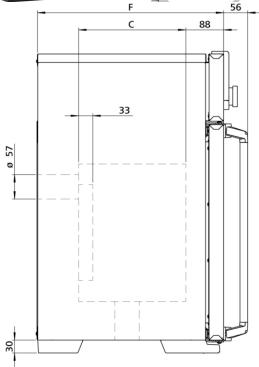

Temperature		
Working temperature range	at least 5 (UN/UNplus/UNm/UNmplus) or 10 (UF/UFplus/UFm/UFmplus) above ambient temperature to +300 °C	
Setting accuracy temperature	up to 99.9 °C: 0.1 / from 100 °C: 0.5	
Setting temperature range	+20 to +300 °C	
Temperature sensor	1 Pt100 sensor DIN class A in 4-wire-circuit	
Control technology		
ControlCOCKPIT	SingleDISPLAY. Adaptive multifunctional digital PID-microprocessor controller with high-definition TFT-colour display	
Timer	Digital backwards counter with target time setting, adjustable from 1 minute to 99 days	
Function SetpointWAIT	the process time does not start until the set temperature is reached	
Calibration	three freely selectable temperature values	
adjustable parameters	temperature (Celsius or Fahrenheit), fan speed, air flap position, programme time, time zones, summertime/wintertime	
Ventilation		
Fan	forced air circulation by quiet air turbine, adjustable in 10 % steps	
Fresh air	Admixture of pre-heated fresh air by electronically adjustable air flap	
Vent	vent connection with restrictor flap	
Communication		
Documentation	programme stored in case of power failure	
Programming	AtmoCONTROL software for reading out, managing and organising the data logger via Ethernet interface (temporary trial version can be downloaded). USB stick with AtmoCONTROL software available as accessory (on demand).	
Safety		
Temperature control	adjustable electronic overtemperature monitor and mechanical temperature limiter TB, protection class 1 according to DIN 12880 to switch off the heating approx. 20°C above nominal temperature	
Autodiagnostic system	for fault analysis	
	101 fault difailysis	
Standard equipment	Tot fault difailysis	
Standard equipment Works calibration certificate	Calibration at +160°C	

Stainless steel interior

Dimensions	$w_{(A)} \times h_{(B)} \times d_{(C)}$: 400 x 320 x 250 mm (d less 39 mm for fan)	
Interior	easy-to-clean interior,made of stainless steel, reinforced by deep drawn ribbing with integrated and protected large-area heating on four sides	
Volume	32	
Max. number of internals	3	
Max. loading of chamber	60 kg	
Max. loading per internal	20 kg	

Textured stainless steel casing

Dimensions	w _(D) x h _(E) x d _(F) : 585 x 704 x 434 mm (d +56mm door handle)	
Housing	rear zinc-plated steel	

Electrical data

Voltage Electrical load	230 V, 50/60 Hz approx. 1600 W	
Voltage Electrical load	115 V, 50/60 Hz approx. 1600 W	

Ambient conditions

Set Up	The distance between the wall and the rear of the appliance must be at least 15 cm. The clearance from the ceiling must not be less than 20 cm and the side clearance from walls or nearby appliances must not be less than 5 cm.
Altitude of installation	max. 2,000 m above sea level
Ambient temperature	+5 °C to +40 °C
Humidity rh	max. 80 %, non-condensing
Overvoltage category	II
Pollution degree	2

Packing/shipping data

Transport information	The appliances must be transported upright
Customs tariff number	8419 8998
Country of origin	Federal Republic of Germany
WEEE-RegNo.	DE 66812464
Dimensions approx incl. carton	w x h x d: 660 x 890 x 650 mm
Net weight	approx. 45 kg
Gross weight carton	approx. 61 kg

Standard units are safety-approved and bear the test marks







