


### IVF module for ICO105med

Ideally suited for in vitro fertilisation.  
Controlled CO<sub>2</sub> and O<sub>2</sub> concentration as  
well as gentle cultivation in Petri dishes in  
separate slide-in units.



On this page, you can find all the essential technical data on the Memmert CO<sub>2</sub> incubator ICOfmed with IVF module. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at [sales@memmert.com](mailto:sales@memmert.com).

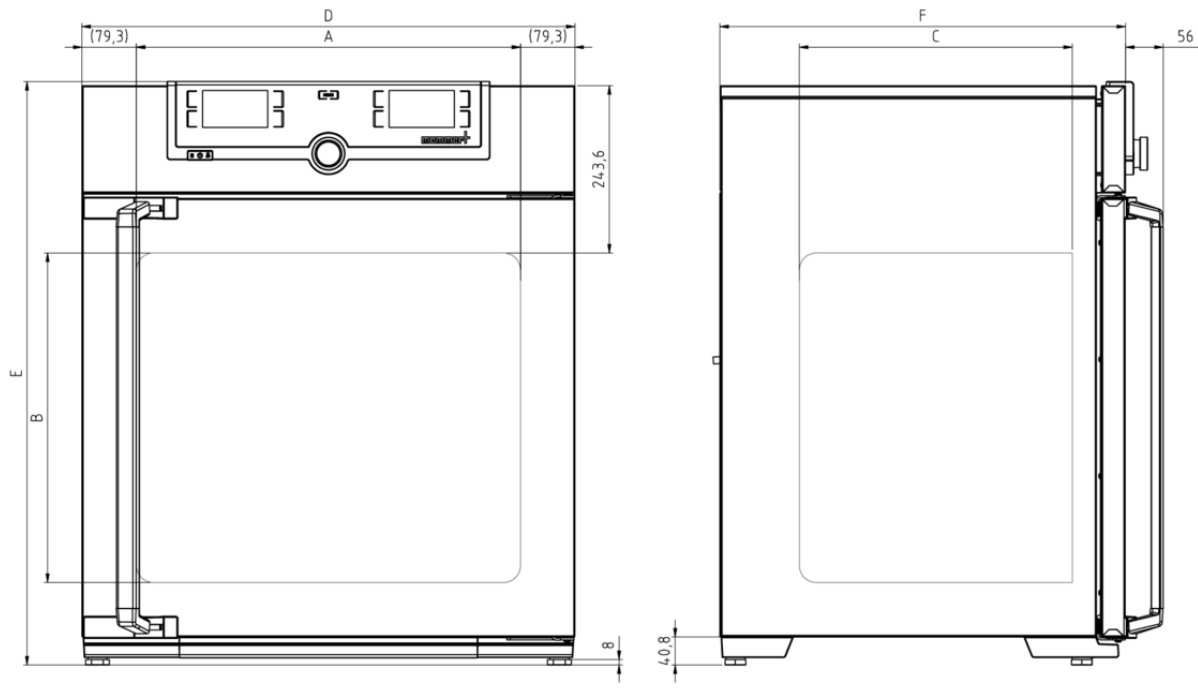

## Temperature

<b>Working-temperature range</b>	5 °C above ambient temperature up to +50 °C Standard sterilisation programme: 60 minutes at 180°C (without removing the sensors)
<b>Setting temperature range</b>	+18 to +50 °C
<b>Setting accuracy temperature</b>	0.1 °C
<b>Temperature sensor</b>	2 Pt100 sensors DIN Class A in 4-wire-circuit for mutual monitoring, taking over functions in case of an error
<b>Temperature variation in chamber</b>	at + 37 °C +/- 0.3 K
<b>Temperature fluctuations with time</b>	at 37 °C +/- 0.1 K

## Humidity

<b>Humidity control (standard)</b>	Humidity limitation thanks to a Peltier element; when water dish is full and inserted, the Peltier element limits the value of relative humidity in the interior to 93 % rh +/- 2.5 %
<b>Setting accuracy humidity</b>	0.5 % rh
<b>Setting range active humidity control (with option K7)</b>	40 to 97 % rh and rh-Off

## Control of standard components

<b>CO<sub>2</sub> control</b>	Digital electronic CO <sub>2</sub> control with dual beam NDIR system, with auto-diagnostic system and acoustic fault indication, barometric pressure compensation
<b>Adjustment range CO<sub>2</sub></b>	0 to 20 % CO <sub>2</sub>
<b>Setting accuracy CO<sub>2</sub></b>	0,1%
<b>Variation in time CO<sub>2</sub></b>	+/- 0.2 % CO <sub>2</sub>
<b>Adjustment range O<sub>2</sub></b>	1 to 20 % O <sub>2</sub>
<b>Setting accuracy O<sub>2</sub></b>	0.1 % O <sub>2</sub>

## Control technology

<b>ControlCOCKPIT</b>	TwinDISPLAY. Adaptive multifunctional digital PID-microprocessor controller with 2 high-definition TFT-colour displays.
<b>Language setting</b>	German, English, Spanish, French, Polish, Czech, Hungarian
<b>Function SetpointWAIT</b>	the process time does not start until the set temperature is reached
<b>Adjustable parameters</b>	temperature (Celsius or Fahrenheit), CO <sub>2</sub> , programme time, time zones, summertime/wintertime

## Communication

<b>Interface</b>	Ethernet LAN, USB
<b>Documentation</b>	programme stored in case of power failure
<b>Programming</b>	AtmoCONTROL software on a USB stick for programming, managing and transferring programmes via Ethernet interface or USB port

## Safety

<b>AutoSAFETY</b>	additionally integrated over- and undertemperature protection "ASF", automatically following the setpoint value at a preset tolerance range, alarm in case of over- or undertemperature, heating function is switched off in case of overtemperature, cooling function in case of undertemperature
<b>Autodiagnostic system</b>	integral fault diagnostics for temperature and CO <sub>2</sub>
<b>Alarm</b>	visual and acoustic

## Heating concept

<b>6 sides</b>	large-area multi-function heating system on four sides with additional door and back heating to avoid condensation
----------------	--

## Standard equipment

<b>Works calibration certificate</b>	Works calibration certificate at +37 °C, 5 % CO <sub>2</sub> and 90 % rh (requires option K7); standard equipment for units with active humidity control
<b>Works calibration certificate</b>	Works calibration certificate at +37 °C, 5 % CO <sub>2</sub> , 90 % rh and 10 % O <sub>2</sub> (requires option K7 and option T6); standard equipment for units with O <sub>2</sub> control
<b>works calibration certificate</b>	works calibration certificate (measuring point chamber centre) at +37 °C, 5 %, 6 % and 7 % CO <sub>2</sub> as well as 90 % rh; 5 % O <sub>2</sub> for IVF unit equipped with option T6
<b>Door</b>	fully insulated stainless steel door with 2-point locking (compression door lock)
<b>Door</b>	inner glass door with opening (8 mm Ø) to take gas sample
<b>Standard equipment</b>	CO <sub>2</sub> connection set: hose with coupling and clamp
<b>Standard accessories</b>	Membrane filter (in order to remove impurities and pollutants, all incoming gases pass through a membrane filter before they reach the chamber)
<b>IVF-module</b>	patented, consisting of 8 slide-in units, a total of 16 special racks with indentations for 16 Petri dishes (60 mm diam.) resp. 32 Petri dishes (35 mm diam.), 2 racks with indentations for 3 special media tubes each; racks with indentations for 4-well dishes on demand;
<b>Humidity control (standart)</b>	Active microprocessor control for humidifying and dehumidifying (40 – 97 % rh), incl. digital indication and auto-diagnostic system ensures even more rapid reaching of set humidity and very short recovery times while avoiding condensate formation. Humidity supply with distilled water (from an external tank) by a self-priming pump; integral bacteria block by generating hotsteam, dehumidifying via sterile filter

### Stainless steel interior

Dimensions	$w_{(A)} \times h_{(B)} \times d_{(C)}$ : 560 x 480 x 400 mm (d less 35 mm for fan)
Interior	material 1.4301 (ASTM 304), corrosion resistant
Volume	107 l
Max. number of internals	6
Max. loading of chamber	90 kg
Max. loading per internal	15 kg

### Textured stainless steel casing

Dimensions	$w_{(D)} \times h_{(E)} \times d_{(F)}$ : 719 x 850 x 591 mm (d +56mm door handle)
Housing	rear zinc-plated steel

### Electrical data

Voltage	230 V, 50/60 Hz
Electrical load	approx. 1300 W
Voltage	115 V, 50/60 Hz
Electrical load	approx. 1300 W

### Ambient conditions

Set Up	The distance between the wall and the rear of the appliance must be at least 15 cm. The clearance from the ceiling must not be less than 20 cm and the side clearance from walls or nearby appliances must not be less than 5 cm.
Ambient temperature	10 °C to 35 °C
Humidity rh	max. 70 %, non-condensing
Altitude of installation	max. 2,000 m above sea level
Overvoltage category	II
Pollution degree	2

### Packing/shipping data

Transport information	The appliances must be transported upright
Customs tariff number	8419 8998
Country of origin	Federal Republic of Germany
WEEE-Reg.-No.	DE 66812464
Dimensions approx incl. carton	w x h x d: 800 x 1030 x 800 mm
Net weight	approx. 75 kg
Gross weight carton	approx. 100 kg

**Standard units are safety-approved and bear the test marks**

